

ALK tel (2p23.2) Orange

FISH Probe
903-7257-102517

BIOCARE
M E D I C A L

Catalog Number: PFA7257V

Description: ALK tel (2p23.2) Orange FISH Probe
Volume: 25 µL


Intended Use:

Analyte Specific Reagent. Analytical and performance characteristics are not established.

Reagents Provided:

The ALK tel (2p23.2) Orange FISH Probe is a fluorescently labeled probe designed to hybridize to ~ 490 kb of the 2p23.2 region of chromosome 2.

ALK tel (2p23.2) Orange probe (60 ng/µl) is provided in FISH Hybridization Buffer (FRR7311A) containing dextran sulfate and formamide.

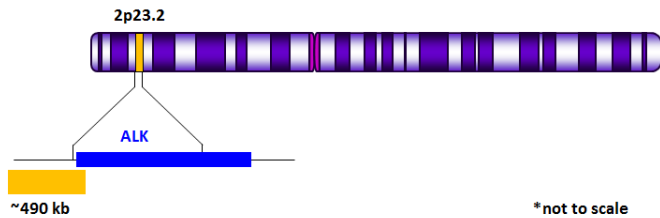

Storage and Stability:

Store probe at -20°C and away from light. The product is stable to the expiration date printed on the label, when stored under these conditions. Do not use after expiration date.

Technical Note:

Biocare Medical FISH probes are optimized to provide the best signal performance using optical filters that can accommodate the excitation/emission wavelengths specified below. Using filters outside these spectral specifications may produce sub-optimal results.

Fluorophore	Excitation (nm)	Emission (nm)
AQUA	432	472
GREEN	498	521
ORANGE	546	575
RED	593	618

Analyte Specific Reagent Note:

1. This product is an Analyte Specific Reagent (ASR). Analytical and performance characteristics are not established.
2. The ALK tel (2p23.2) Orange FISH probe may cross-hybridize to other areas on the human genome. It is the responsibility of the user to validate any test for its specific use.

Precautions:

1. This product contains formamide, which may be toxic. Formamide may cause serious eye damage or reproductive toxicity. It may also cause irritation by inhalation or skin contact. Avoid any direct contact exposure to reagent. Take appropriate protective measures (use disposable gloves, protective glasses, and lab garments).
2. Specimens, before and after fixation, and all materials exposed to them should be handled as if capable of transmitting infection and disposed of with proper precautions. Never pipette reagents by mouth and avoid contacting the skin and mucous membranes with reagents and specimens. If reagents or specimens come in contact with sensitive areas, wash with copious amounts of water¹.
3. The SDS is available upon request and is located at <http://biocare.net>.

References:

1. Clinical and Laboratory Standards Institute (CLSI). Protection of Laboratory workers from occupationally Acquired Infections; Approved Guideline-Fourth Edition CLSI document M29-A4 Wayne, PA 2014.

 Biocare Medical

60 Berry Drive
Pacheco, CA 94553

USA

Rev. 062117

Tel: 800-799-9499 | www.biocare.net | Fax: 925-603-8080